

QCT PRO™


PACS Option

FEATURES

- Receipt via DICOM of data from multiple CT scanners
- Export of QCT PRO Spine and/or Hip BMD reports to DICOM destinations
- Multiple methods for selecting the destination(s) for a report, including:
 - static list of one or more destinations,
 - selection from a drop-down list,
 - automatic delivery based on CT scanner ID associated with an analyzed data set

The QCT PRO PACS option provides integration of QCT PRO into your local PACS solution. Integration allows you to exploit your local PACS infrastructure to work more efficiently by providing a paperless report that can be electronically archived and retrieved, reviewed locally or remotely, delivered in electronic form to your patients or physicians, or handled in other manners supported by your PACS.

After installation, PACS export is automatic and occurs during report printing.

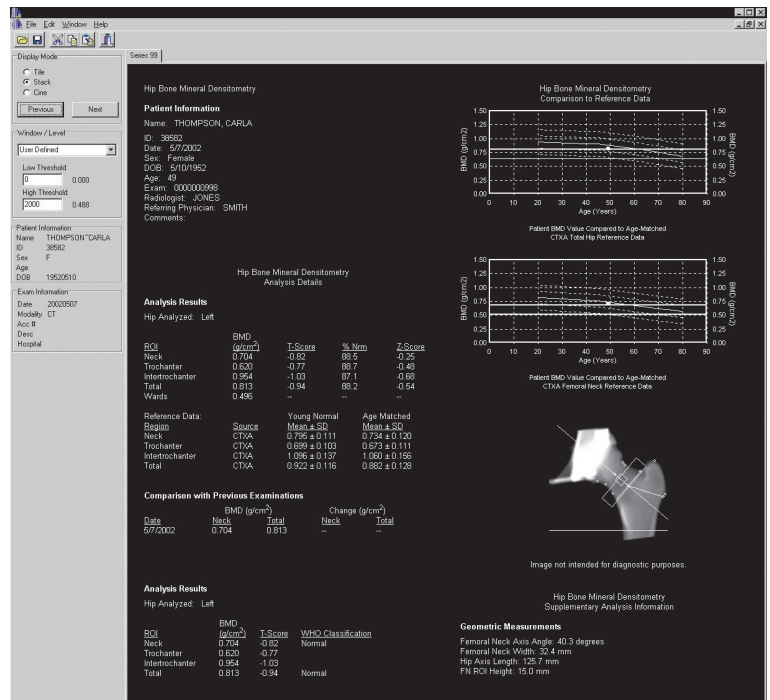

Exact display of information depends on your PACS system.

Mindways Software, Inc.
3001 South Lamar Blvd.
Austin, TX 78704 USA

Tel: 512 912 0871
Email: info@qct.com
www.qct.com