




Extract BMD measurements from routine\* CT or PET/CT abdomen studies without extra scans or exposure

# SlicePick™

\*No I.V. contrast

SlicePick is a QCT PRO™ BMD image extraction utility. Simply include a QCT calibration phantom in the scan and use SlicePick to extract the images for BMD analysis with QCT PRO. The phantom does not interfere with the abdomen study.

SlicePick generates AP and lateral projection images from volumetric CT data sets. The vertebral bodies to be measured are enclosed within a box and extracted, ready for analysis with QCT PRO. QCT PRO's 3D volumetric technology enables precise ROI placement without the gantry angulation required with conventional 2D QCT.

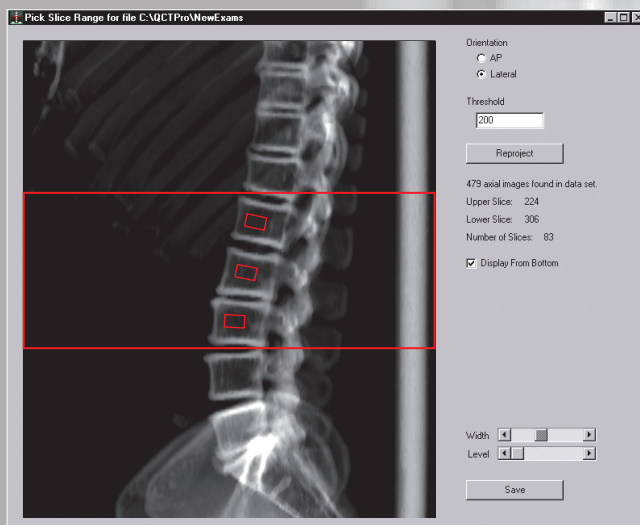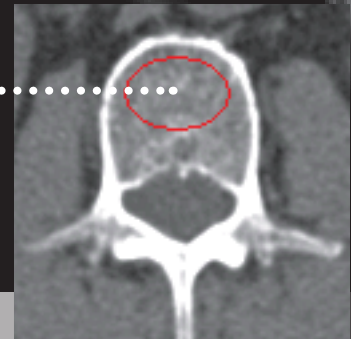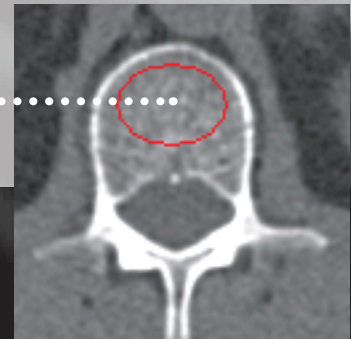

Mindways Software, Inc.  
3001 S Lamar Blvd Ste 302  
Austin TX 78704-4799

Tel: 512 912 0871  
Email: [info@qct.com](mailto:info@qct.com)  
[www.qct.com](http://www.qct.com)